

ACCENT® 980

Laser Lipolysis Accentuation

This 980 nm diode laser is the perfect wavelength for minimally invasive laser assisted lipolysis. Energy and heat are transmitted through an optical fiber and delivered to the subcutaneous fat for ease of removal and simultaneous collagen stimulation



Accent 980 Features


- FDA approved for laser assisted lipolysis
- Procedure done under local anesthetic
- Collagen stimulation for improved skin laxity
- Faster and less traumatic aspiration process
- Less downtime and bruising
- Compact carry case with impressive portability
- Popular patient procedure
- Procedure produces high revenue

Accent 980 Benefits

- No off label advertising
- Increases safety and greater patient appeal
- Tuck and lift is less necessary after procedure
- Requiring less physician exertion
- Quick recovery (up and around in 24-48 hrs)
- Easily transportable between multiple offices
- Won't lose patients to competition
- Increase total practice revenue with a high ROI



"The Accent 980 adds a whole new dimension in the way we treat increased volume and skin laxity. Additionally, compared to traditional liposuction, which involves seven to ten days of downtime, the Accent 980 is minimally invasive with minimal pain or downtime and 90% less bruising. Even in this challenging economic environment, these types of treatments are big sellers. Results are also tremendous. It's a home run." Dr. Karen Harkaway, Board Certified Dermatologist from Delran, NJ

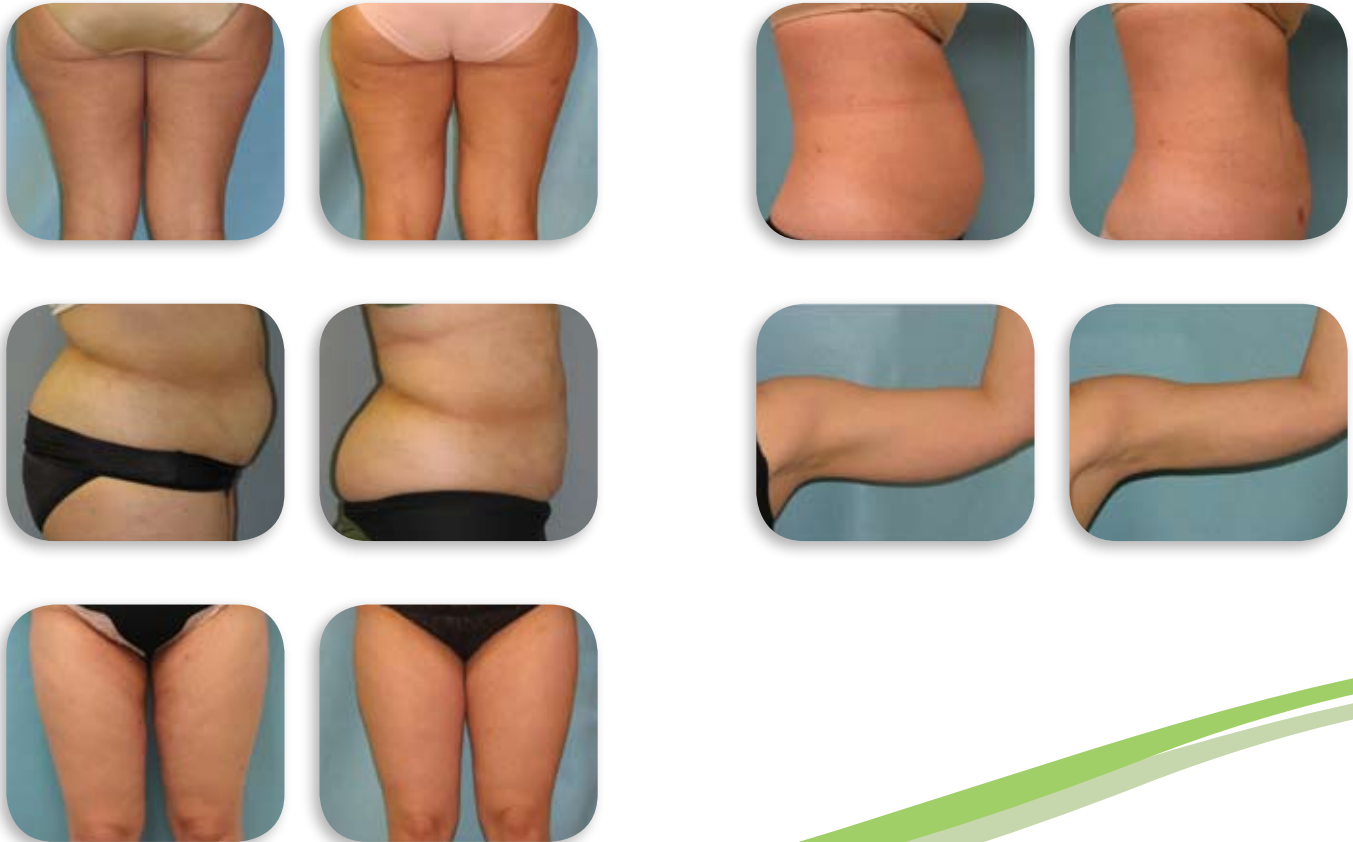

Courtesy of Karen Harkaway, MD

Specifications

Wavelength	980 nm
Power:	10 watts
Pulse length:	1 ms to cw
Pulse interval:	2 ms to 1 sec
Fiber size:	400 microns
Aiming beam:	640 nm, 2 mW
Dimensions:	12 x 21 x 10 cm
Weight:	10 Lbs
Rechargeable battery or 115 AC	

©2008 Alma Lasers, Ltd. All rights reserved. Alma Lasers Ltd., its logo, and Pixel® are the trademarks of Alma Lasers, Ltd. Product specifications are subject to change without notice.

Alma Lasers®
Wellbeing Through Technology®

US Headquarters
485 Half Day Road
Suite 100
Buffalo Grove, IL 60089
Tel 866-414-ALMA
Fax 224-377-2050
contact@almalasers.com

International Headquarters
14 Halamish st.
Caesarea Industrial Park
Caesarea, 38900 Israel
Tel +972-4-627-5357
Fax +972-4-627-5368
info@almalasers.com